

TRAINING UPDATE

Lab Location:	GEC, SGAH & WAH	Date Distributed:	4/15/2014
Department:	Mgmt, QA & LIS teams, Technical Specialist	Due Date:	4/30/2014
		Implementation:	5/1/2014

DESCRIPTION OF PROCEDURE REVISION

Name of procedure:
Process and Equipment Validation Protocol GEC / SGAH / WAH.QA46 v1
New form added: Checklist for New/Changed Procedure Implementation (AG.F279)
Description of change(s):
Section 5: add process for using Checklist for New/Changed Procedure
This revised SOP will be implemented on May 1, 2014

Document your compliance with this training update by taking the quiz in the MTS system.

Approved draft for training all sites (version 1)

Non-Technical SOP

Title	Process and Equipment Validation Protocol	
Prepared by	Leslie Barrett	Date: 7/27/2012
Owner	Cynthia Bowman-Gholston	Date: 7/27/2012

Laboratory Approval		
Print Name and Title	Signature	Date
<i>Refer to the electronic signature page for approval and approval dates.</i>		
Local Issue Date:		Local Effective Date:

Review:		
Print Name	Signature	Date

TABLE OF CONTENTS

1. PURPOSE.....	3
2. SCOPE	3
3. RESPONSIBILITY.....	3
4. DEFINITIONS.....	3
5. PROCEDURE.....	4
6. RELATED DOCUMENTS	5
7. REFERENCES	5
8. REVISION HISTORY.....	5
9. ADDENDA AND APPENDICES	5

1. PURPOSE

This procedure describes the process of validation and revalidation of new and revised processes or procedures, and new or repaired equipment.

2. SCOPE

This procedure applies to all non-computer processes or procedures and equipment that require validation.

3. RESPONSIBILITY

- **Testing Personnel** will:
 - Follow the Validation Protocol as designed.
 - Document all steps of the validation protocol.
- The **Operations Director, Manager, or Supervisor** (or designee) will:
 - Ensure compliance with this procedure in his/her department.
 - Document all steps of the validation protocol.
 - Ensure that all staff is appropriately trained.
- The **Quality Assurance Department** will:
 - Manage validation documents on the electronic document control system
- The **Medical Director** will:
 - Approve the initial document and any subsequent revisions.
 - Approve all validations prior to implementation.
 - Approve validation protocols.

4. DEFINITIONS

Installation Qualification (IQ): ensures equipment is installed per requirements

- Equipment is installed per manufacturer's specifications
- Equipment operates within the limits required for the process and per manufacturer specifications

- Support utilities (water, electric, air, temperature) must meet manufacturer's specifications

Operational Qualification (OQ): ensures equipment can do what it is supposed to do. Examples: test alarms, check temperatures of empty refrigerator, test RPM/RCF and timers of centrifuge, QC results, etc.

Performance Qualification (PQ): ensures equipment does perform as expected. Test using normal samples under normal and extreme conditions. Must be performed by staff who will utilize the equipment.

Examples: fill the refrigerator and ensure it still maintains temperature in all areas of the refrigerator; ensure filling the refrigerator does not inhibit air flow; fill a centrifuge and make sure it runs appropriately with maximum and/or minimum load, etc.

Validation: “verification, where the specified requirements are adequate for intended use.” (reference: “International vocabulary of metrology, Basic and General concepts and associated terms (VIM) 2008”).

Verification: “provision of objective evidence that a given item fulfills specified requirements.”

5. PROCEDURE

1. The Supervisor or designee develops, writes, and submits validation protocol to include laboratory regulations, accreditation standards, and manufacturer’s instructions. Refer to addenda B ~~and C~~ for ~~form and~~ guidelines.

Notes:

- [An Equipment and Process Validation Protocol Form must be used to document the plan. Refer to Related Documents](#)
 - For validation of test methods refer to the procedure Laboratory Method Validation for Quantitative and Semi-Quantitative Methods, QDQC710.
2. The Laboratory Director reviews and approves draft protocol, or revises as needed.
 3. Installation Qualification and Operational Qualification are performed and documented by the appropriate personnel. This may include Biomedical personnel, manufacturer’s representative, in-house or other designated staff.
 4. The Supervisor reviews and approves IQ and OQ documentation. Corrective action is documented as necessary.
 5. The Supervisor or designee writes draft procedure(s) and training document. Pre-validation training is performed and documented for appropriate staff.
 6. Assigned validation staff performs the Performance Qualification and validation. Data is collected and documented on appropriate forms (ie, QC, Maintenance logs)

7. The Checklist for New/Changed Procedure Implementation is used to document completion of key tasks for technical processes (example: LIS changes, QC definitions, communication of changes to reference ranges or critical values, etc).
 - The form prompts for dates and initials for each task. Refer to Related Documents
 - The project champion is responsible for assuring each applicable item is completed, including appropriate documentation.
 - The form must be submitted with the final validation documentation when changes have the potential to impact analytical tests or test systems.
8. The Technical Supervisor and Medical Director review validation data and approve/disapprove validation via electronic document control system. The procedure(s) is finalized and approved via electronic document control system.
9. The Supervisor or assigned validation staff train remaining staff and document.
10. The process or equipment is implemented.

6. RELATED DOCUMENTS

- Laboratory Method Validation for Quantitative and Semi-Quantitative Methods, QDQC710
- Process/Equipment Validation Protocol, Transfusion Service, QDHBB303
- SmartSolve® (Pilgrim) EDCS: Managing New, Revised, Expire and Recurring Review of Documents, QA procedure
- Equipment and Process Validation Protocol Form (AG.F204)
- Checklist for New/Changed Procedure Implementation (AG.F279)

7. REFERENCES

Process/Equipment Validation Protocol, Transfusion Service, QDHBB303

8. REVISION HISTORY

Version	Date	Reason for Revision	Revised By	Approved By
000	3/25/14	Section 5: add process for using Checklist for New/Changed Procedure Section 6: update EDCS title, add Checklist form Section 9: form moved to section 6 Footer: version # leading zero's dropped due to new EDCS in use as of 10/7/13.	L Barrett	C Bowman


9. ADDENDA AND APPENDICES

- A. Process Validation Flow Chart
- B. Validation Protocol Format
- ~~C. Equipment and Process Validation Protocol Form (see Attachment Tab of Infocard)~~

Form revised 3/31/00

ADDENDUM A

PROCESS VALIDATION FLOW CHART

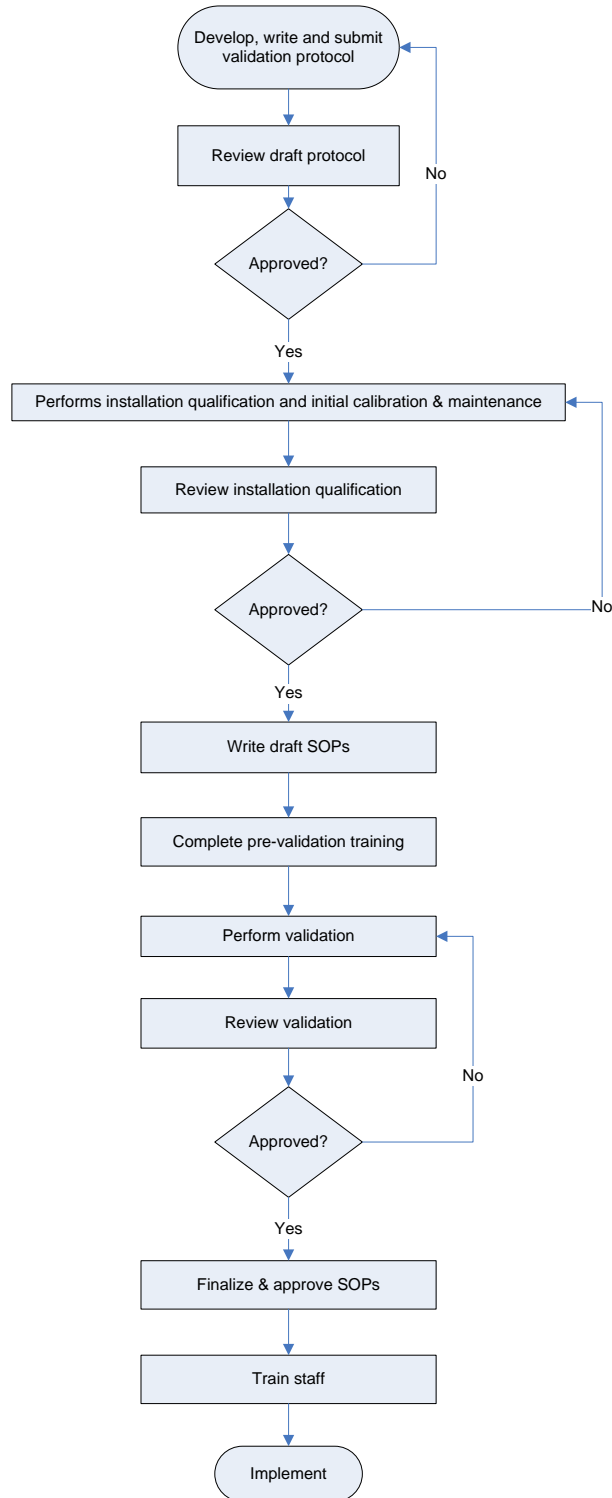

Form revised 3/31/00

ADDENDUM B VALIDATION PROTOCOL FORMAT

The following items are provided as guidelines for a validation protocol.

A. Title and protocol approvals

This section contains the protocol title and a space for a signature approving the protocol.

B. Purpose

1. Identify what equipment or process is being addressed in the validation.
2. Identify the reason the validation or revalidation is required

C. Requirements and Specifications

1. Identify the performance requirements for the equipment or process with the following:
 - a. Equipment capabilities and parameters
 - b. Regulatory requirements and industry standards associated with the equipment or process
2. References, (operator's manual, package insert, etc)

D. Test Plan

1. Installation Qualifications
 - a. Identify all steps required to verify that the equipment was properly installed in environmental conditions that meet the manufacturer's specification.
 - b. Provision for sign-off that the installation qualification has been met.
2. Operational Qualifications
 - a. Identify all steps required to verify that the equipment or process operates within established limits and specifications supplied by the manufacturer.
3. Performance Qualifications
 - a. Identify all steps required to verify that the equipment performs as expected for its intended use in the processes established by the facility and that the output meets the facility's specification.
4. General requirements
 - Each test should compare the expected results from the requirements against the actual results.
 - All steps should provide a space for the person performing the validation to initial and date.
 - The validation plan should specify the number of process runs, when applicable, to provide an accurate measure of variability among successive runs.

Checklist for New/Changed Procedure Implementation

Test/ Instrument:

Checklist Champion:

(Champion oversees checklist is followed and obtains all necessary documentation and signatures prior to implementation.)

Est. start date:

Est. Completion:

Description:

Checklist		
Proposed change approved	Initial	Date
Technical Director	<input type="text"/>	<input type="text"/>
Approved Valication Plan	Initial	Date
(Linearity,AMR,CRR,Carryover,etc)		
Technical Director	<input type="text"/>	<input type="text"/>
Medical Director	<input type="text"/>	<input type="text"/>
Establishment and verification of method performance specifications	Initial	Date
Analytical sensitivity	<input type="text"/>	<input type="text"/>
Precision	<input type="text"/>	<input type="text"/>
Accuracy	<input type="text"/>	<input type="text"/>
Approved SOP	Initial	Date
Procedure approved and under document control	<input type="text"/>	<input type="text"/>
Forms approved and under document control	<input type="text"/>	<input type="text"/>
Test equipment,reagents, and supplies on hand	Initial	Date
	<input type="text"/>	<input type="text"/>
Equipment maintenance and function checks	Initial	Date
	<input type="text"/>	<input type="text"/>
Proficiency Testing	Initial	Date
Information provide to QA _____	<input type="text"/>	<input type="text"/>
Additional actions required (specify)	Initial	Date
	<input type="text"/>	<input type="text"/>
IT	Initial	Date
LIS change request submitted,must include the following:	<input type="text"/>	<input type="text"/>
CDM,CPT,LOINC	<input type="text"/>	<input type="text"/>
Report Format	<input type="text"/>	<input type="text"/>
AHC Cerner changes	<input type="text"/>	<input type="text"/>
Verify reference Ranges are identical in procedure,instrument	<input type="text"/>	<input type="text"/>
Validate instrument interface	<input type="text"/>	<input type="text"/>
Validate Computer Interfaces-LIS,HIS,POCT	<input type="text"/>	<input type="text"/>