



Atrium Health

Wake Forest Baptist

This material is for demonstration purposes only and does not contain real patient information.




WakeOne Winter Update 2021

BEAKER

Beaker – Anatomic Pathology

These updates are specific to Anatomic Pathology.

Edit Results Before Specimens Are Marked Gross Done

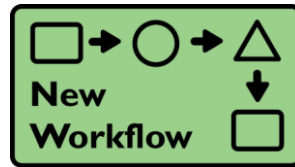

- Users entering or editing results for a case no longer needed to wait for grossing assistants to mark a specimen done before they edit case information in the Result Editor pane of the Outstanding List. They can see the in-progress Gross Description for specimens and edit other fields in the meantime.
- As before, specimens that no users have started grossing do not appear.

The screenshot shows the 'Tissue Exam' software interface. On the left, there is a 'Clinical Information [1]' section with a text box containing 'Sections through the liver show distorted hepatic architecture.' Below it is a 'Gross Description [2]' section with two entries: 'A. Liver.' and 'B. Liver.'. Entry 'A' is highlighted with a red box and contains the text: 'The specimen is received fresh. It consists of a total hepatectomy, 380 grams and 12.5 x 10.0 x 7 cm. The capsular surface is smooth and pink.' On the right, there is a 'Specimens' table with columns for ID, Protocol, Source, Description, Collected, and Received. The table contains three rows: 'A Liver, Biopsy', 'B Liver, Biopsy', and 'C Liver, Biopsy'. The 'Received' column for row 'A' has a red box around a clock icon, while row 'B' has a green checkmark. Below the table is a 'Questions' section with a table for 'Collection Question' and 'Specimen Source Modifier (free text)'. The top navigation bar includes 'Snapshot', 'Forms', 'Coding', and 'QA'. The left sidebar has 'Routine' and 'Resulted' buttons.

ID	Protocol	Source	Description	Collected	Received
A	Liver, Biopsy	Liver		02/19/21 0837	02/19/21 0837
B	Liver, Biopsy	Liver		03/01/21 0704	03/01/21 0704
C	Liver, Biopsy	Liver		03/01/21 0704	03/01/21 0707




New Specimen Numbering Scheme


- Beaker can use placeholder specimens to reserve an ID when specimens arrive in the lab out of sequence so that the ID of the specimen collected in the OR and the ID of the specimen in Case Builder match more often. A specimen is considered a placeholder specimen until the point at which it is acknowledged as part of a case in Case Builder.
- When a case has more than 26 specimens, **the ID sequence now goes from Z to AA, AB instead of AA, BB as it did previously** so that specimen IDs match the naming convention of OR specimens. This change means that when additional specimens are added to a case that already had some double-letter labeled specimens, the new specimens fill into the previous list alphabetically instead of being appended.

5460740-Y
5460740-Z
5460740-AA
5460740-AB
5460740-AC
5460740-AD

Jump to Case QA in the Outstanding List

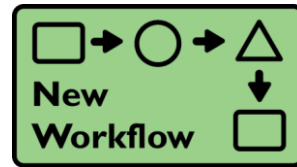

- When you click the Case QA From the following places, you are now redirected to the Case QA tab in the right sidebar of the web-migrated Outstanding List instead of to the separate activity so that your workflow is consistent no matter where you access Case QA from.
 - Case Inquiry
 - Reporting Workbench
 - Follow-up Work List

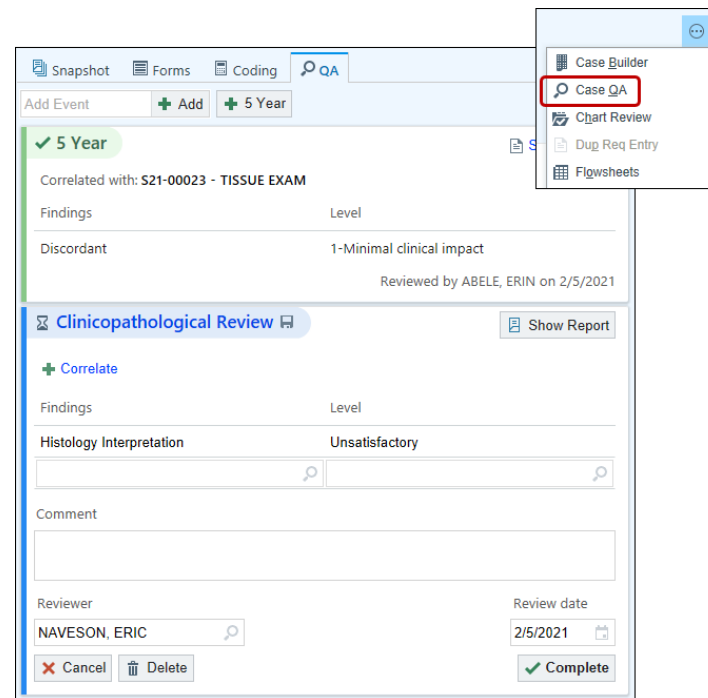

Beaker – Clinical Pathology

These updates are specific to Clinical Pathology.

The Web-Migrated Outstanding List is Coming to Clinical Pathology



- The Outstanding List for Clinical Pathology will be changing to the Web-Migrated Outstanding List.

Outstanding List

M	ID	Container ID	L... Name	F P	Test	A D L R H Q..	TAT Remain...	Status
	21W-299CH0003	9000152406	Fake-Patient, Jessica Beaker		T4		1h 53m	On Test
	21W-299CH0003	9000152406	Fake-Patient, Jessica Beaker		Free T4		1h 53m	On Test
	21W-299CH0003	9000152406	Fake-Patient, Jessica Beaker		TSH, 3rd GENERATION		1h 51m	On Test
	21W-299CH0002	9000152405	Fake-Patient, Jessica Beaker		Complete Metabolic Panel		48m	On Test
✓	21W-299CH0001	9000152401	Fake-Patient, Jessica Beaker		Potassium		1h 40m	On Test
	21W-299CH0001	9000152401	Fake-Patient, Jessica Beaker		Magnesium		1h 40m	On Test
	21W-299CH0001	9000152401	Fake-Patient, Jessica Beaker		FREE T3		5h 40m	On Test
	21W-295CH0001	9000152369	Grandcentral, Corona		Kidneyintex		-3d 12h 26m	On Test
	21W-294CH0001	9000152284	Grandcentral, Arizona		Kidneyintex		-4d 15h 6m	On Test
	21W-293CH0006	9000152171	Grandcentral, Wyoming		Kidneyintex		-5d 16h 34m	On Test
	21W-293CH0005	9000152170			intex		-5d 16h 34m	On Test
	21W-293CH0001	9000152165			intex		-5d 18h 57m	On Test
	21W-291CH0008	9000152028			nel		-7d 17h 52m	On Test
	21W-252CH0003	9000148981			Urine		-46d 17h 17m	On Test
	21W-252CH0002	9000148980	Test, Emily		Urea Nitrogen, Urine		-46d 17h 17m	On Test
	21W-252CH0002	9000148980	Test, Emily		Protein, Urine		-46d 17h 17m	On Test
	21W-252CH0002	9000148980	Test, Emily		Phosphorus, Urine		-46d 17h 17m	On Test
	21W-252CH0002	9000148980	Test, Emily		Sodium/Potassium, Urine		-46d 17h 17m	On Test
	21W-252CH0002	9000148980	Test, Emily		Creatinine, Urine		-46d 17h 17m	On Test
	21W-252CH0002	9000148980	Test, Emily		Citrate, Urine		-46d 17h 17m	On Test
	21W-252CH0002	9000148980	Test, Emily		Chloride, Urine		-46d 17h 17m	On Test
	21W-252CH0002	9000148980	Test, Emily		Calcium, Urine		-46d 17h 17m	On Test
	21W-252CH0002	9000148980	Test, Emily		Uric Acid, Urine		-46d 17h 17m	On Test
	21W-252CH0002	9000148980	Test, Emily		Albumin, Urine		-46d 17h 17m	On Test
	21W-252CH0001	9000148070	Test, Emily		Magnesium Urine		-46d 17h 17m	On Test

Patient Header

Fake-Patient, Jessica... 21W-299CH0001
 MRN 6159997 Blood
 F, 34 years, 10/10/1987, Collected 10/26/2021 0905
 Received 10/26/2021 0905
 Instrument ID 9000152401
 Draw Type Venipuncture

Results Editor

Potassium
 Holds, Flags, FYIs, & Follow Ups
 Hard Draw
 CBC needs to be collected in a Lt. Blue tube, due to platelet clumping.

Res	Component	Value	Units	! Δ L IE R	Ref. Range	Method	Trend	PV
1	Potassium		MMOL/L		3.5-5.3	WCCHR AU5812		2.8

Comment: 1
 1 Hemolysis Serum Index WCCHR AU5812
 1 Icterus Serum Index WCCHR AU5812
 1 Lipemia Serum Index WCCHR AU5812
 Method: WCCHR AU5812 Last received: 10/26/2021 0905

Sidebar

Questions
 Order Question Answer
 Charge To Research Study (To be answered by Research Providers ONLY!)

Previous Results
 8/14/2020 7/10/2020 6/20/2020 6/18/2020
 Component 0924 1030 1445 1437
 Potassium 2.8 5.5 3.2 4.0
 MMOL/L¹ MMOL/L MMOL/L MMOL/L
 Component Comments
¹ CALLED TO AND READBACK BY: Felicity Smoak, MD
 Electronically signed by: Tech Beaker Lab 8/15/2020 1:48 AM

Batch Summary

Test	Batch	Type	Status
Magnesium	<on queue>		
Potassium	<on queue>		

Alternate Specimen IDs

Internal Specimen	74269
21W-299CH0001.1	9000152401


Linked Orders
 None

Expand/Collapse

More Info about the Web-Migrated Outstanding List – 1 of 2



- This updated activity can search the database of tests to generate simple reports and includes a personal work list called Pinned to which users can add tests and cases they intend to work on.
- In the web-migrated Outstanding List, the upper list section is replaced with filter buttons users can click to see only overdue tests and cases. These filters appear by default and require no setup.
- You will need to click the arrow to see the STAT, Overdue and Resulted lists.

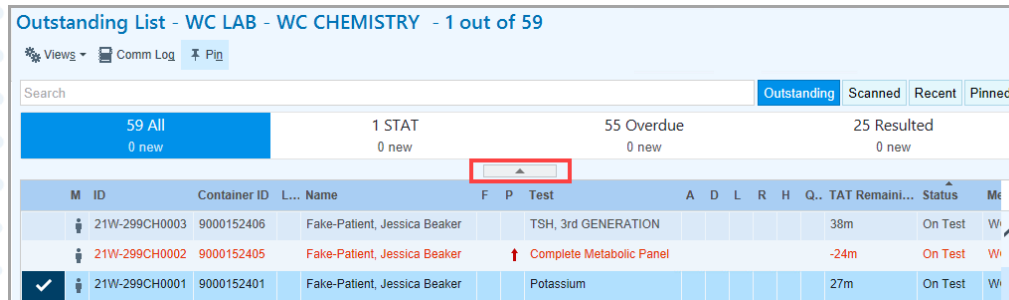

The screenshot shows the 'Outstanding List - WC LAB - WC CHEMISTRY - 1 out of 59' interface. It features a search bar, filter buttons for 'Outstanding', 'Scanned', 'Recent', and 'Pinned', and a summary row with counts: 59 All (0 new), 1 STAT (0 new), 55 Overdue (0 new), and 25 Resulted (0 new). Below is a table with columns for M, ID, Container ID, L... Name, F, P, Test, A, D, L, R, H, Q., TAT Remain..., Status, and Me. A red box highlights an upward-pointing arrow in the 'P' column of the second row.

M	ID	Container ID	L... Name	F	P	Test	A	D	L	R	H	Q.	TAT Remain...	Status	Me
	21W-299CH0003	9000152406	Fake-Patient, Jessica Beaker			TSH, 3rd GENERATION							38m	On Test	W
	21W-299CH0002	9000152405	Fake-Patient, Jessica Beaker		↑	Complete Metabolic Panel							-24m	On Test	W
✓	21W-299CH0001	9000152401	Fake-Patient, Jessica Beaker			Potassium							27m	On Test	W

- When a selected result is verified or otherwise removed from the Outstanding List, no new row in the list is selected and the user sees a message in the central result editing pane that they can select a row to get started. This behavior helps ensure that users are consciously selecting the next test or case they want to work on and are verifying details as needed.



More Info about the Web-Migrated Outstanding List – 2 of 2

- You will still result your tests using the Result Editor. You will access Result Editor by scanning the Specimen ID, double-clicking the patient row or right-click the patient row and select Result Editor.

	21W-299CH0002	9000152405	Fake-Patient, Jessica Beaker	↑ Complete Metabolic Panel	
<input checked="" type="checkbox"/>		21W-299CH0001	9000152401	Fake-Patient, Jessica Beaker	Potassium
	21W-299CH0001	9000152401	Fake-Patient, Jessica Beaker	Magn	
	21W-299CH0001	9000152401	Fake-Patient, Jessica Beaker	FREE	
	21W-295CH0001	9000152369	Grandcentral, Corona	Kidne	
	21W-294CH0001	9000152284	Grandcentral, Arizona	Kidne	
	21W-293CH0006	9000152171	Grandcentral, Wyoming	Kidne	

- Inquiry
- Specimen Update
- Comm Log
- Add Follow-up Task
- Results**
- Pin

- Users have the option to lock the width of the center result entry pane in the web-migrated Outstanding List so that the pane mimics as consistently as possible the way the results will appear on the result report.

Instrument ID
9000152401
Draw Type
Venipuncture

- Remove from Scanned on Verify
- Show Component Hotkeys
- Show Component Rulers
- Lock Resulting Pane Width


- Add Follow-up Task
- Add to Batch
- Add to Packing List
- Additional Settings**
- Cancel
- CC Results
- Chart Review

AP/CP - Cross Lab Updates

These updates will affect users from **both** the AP and CP labs.

More Accessible Reevaluate Result Checking and Report Flags Actions

- The Reevaluate Result Checking and Reevaluate Component Report Flags actions, which were previously located in Result Entry and Verification, are now also available in the web-migrated Outstanding List actions menu, saving you time when you need to perform these tasks.

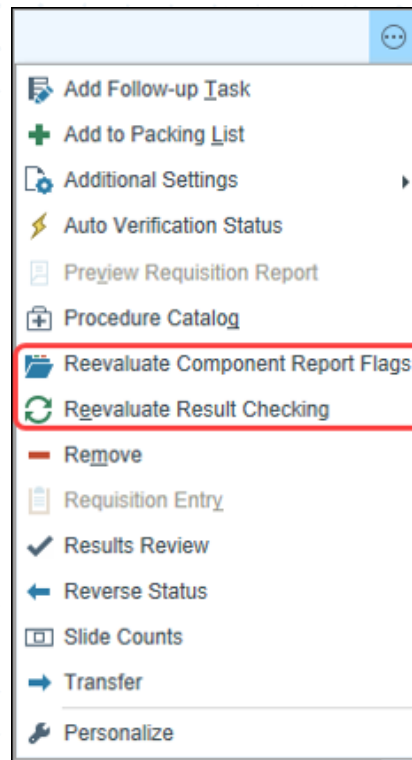



Break the Glass to See Restricted Patients in the Lab Past Due Orders Report

- When you review past due orders that need to be collected, you might notice patient names now appear as "<Restricted>". When needed, click Collect Specimens or the Access patient information link in the Order Summary report to break the glass for the patient so you can see the Order Summary report and document collection. Click Refresh Selected for the selected row to update the patient name column.

Lab Past Due Orders Report [2906437] as of Mon 4/19/2021 12:41 PM

Collect Specimens Saved Views - Original View

Detail List By Ordering Department By Ordering Provider By Procedure By Collection Date

Filter Re-run Report Refresh Selected Select All

Exp Col Date/Time	Order Date	Priority	Patient	Procedure	Lab Order Department	Ordering Provider	Authorizing Provider	Qur Specimen Type	Order Type	Class	Order ID
10/14/2020 0600	10/13/2020	Routine	<Restricted>	PROTIME-INR [LAB320]	EMH MED SURG	Physician Family Medicine	Physician Family Medicine, MD	Blood	Lab	Lab Collect	1065274
10/15/2020 0600	10/14/2020	Routine	Willow, Eric E (201776)	PROTIME-INR [LAB320]	EMH MED SURG	Physician Family Medicine	Physician Family Medicine, MD	Blood	Lab	Lab Collect	1044772
10/15/2020 1439	10/15/2020	Pending Discharge	Beacon, Natalie (203015)	TISSUE EXAM [LAB1126]	EMH HEM ONC	Physician Surgery	Physician Surgery, MD	✓ Tissue	Pathology and Cytology	Unit Collect	1065642

Order Summary

Break-the-Glass

This patient record has been deemed confidential and restricted. [Access patient information](#)

Update Resources

For both AP and CP labs.

WakeOne Update Material

Beaker – Anatomic Pathology



After reviewing this document, use it as a guide to practice the Winter 2021 version features in the **WakeOne Playground**. Playground refreshes nightly and does not affect the live environment.

WakeOne PLY User IDs

Techap##
suppap##
cytoap##
scytoap##
pathap##
refeap##
superap##
techapa##

Password for all Beaker Users: **train**

Use the link below to access the patient information for PLY

[WakeOne Playground Logins - AP](#)

WakeOne Update Material

Beaker – Clinical Pathology



After reviewing this document, use it as a guide to practice the Winter 2021 version features in the **WakeOne Playground**. Playground refreshes nightly and does not affect the live environment.

24 Playground User IDs Available

01 – 24

labphleb##

labasst##

labtech##

Password for all Beaker Users:

train

First Names

Last Names

Hannah

Blushwine

Gabriel

Chenin

Doris

Drywhite

Irene

Fortwine

Juanita

Icewine

Alicia

Merlot

Update Resources

- Learning Home Dashboards
 - Press F1 while logged into WakeOne
- [WakeOne Playground Logins – AP](#)
- [WakeOne Playground Logins – CP](#)

- The WakeOne Playground (PLY) will be updated to the new version – May 2021, and be [available on 11/22/2021](#).